

23OCT2025

URGENT MEDICAL DEVICE CORRECTION
IV ADMINISTRATION SETS & PUMP ADMINISTRATION SETS
(INFUSOMAT® SPACE; OUTLOOK®; & VISTA BASIC)
Occlusion and Backflow

Dear Valued Customer:

The purpose of this letter is to inform you that B. Braun Medical Inc. (BBMI) is issuing a voluntary Urgent Medical Device Correction for gravity IV administration sets as well as pump administration sets utilized with BBMI's Infusomat® Space Large Volume Pump, Outlook Pump, and Vista Basic Pump due to the potential for backflow of medication from secondary (piggyback) IV containers into primary IV containers or other connections and the inability to prime (occlusion).

Affected Products:

Please refer to Attachment 1 for a list of all affected products. The depth of this voluntary correction is to the hospital/healthcare facility level. This correction is being conducted with the knowledge of the United States Food & Drug Administration.

Reason for the Correction:

An investigation of customer complaints has identified that the check valve components of the affected articles may become stuck in the open or closed positions. If the valves remain in the open position, then backflow of medication may occur (e.g. secondary container fluid backs-up into primary container during secondary "piggyback" administration). If the valves remain in the closed position, then difficulty or inability to prime the line may be observed.

Risk to Health:

To date, there have been no reports of serious injury or death associated with this issue. The reported complaint rate is approximately 11 complaints per million units distributed.

- If the user is unable to prime the infusion set, there is a risk of delay in medication delivery while a replacement administration set is obtained.
- If the secondary medication backs up into the primary IV container, there is a risk that the intended dose of medication is either not delivered or is delivered over an unintentionally extended timeframe (coupled with the primary infusion).
 - If the error is identified, additional infusions may be required, with the resultant final dose subject to error given the difficulty quantifying how much of the original infusion was delivered to the patient.
 - If the error is not identified (more likely with clear infusions), underdosing may occur;
 - overdosing is also possible with repeated scheduled infusions when secondary medication that has accumulated in the primary IV container is inadvertently infused into the patient at the primary infusion rate.

Depending on the time delay in infusion, the dosage error, the nature of the infusion, and the medical status of the patient, varying clinical effects may result, up to and including serious injury or death. Note that per the product Instructions for Use (IFU) secondary infusions should not be used for critical infusions.

B | BRAUN

Actions Taken by BBMI:

The reported events have been linked to processing variation from the supplier of the internal disk of the valves. BBMI conducted an on-site audit of the disk supplier and has identified several factors which may influence the performance of the disk. BBMI is currently working with the supplier to implement corrective actions to prevent recurrence of these issues. Until such time as lots with these corrective actions are available, users are asked to follow the "Interim Measures for Users." A follow-up communication will be provided to users when lots with these corrective actions are in place and the interim measures can be stopped.

Interim Measures for Users:

1. Secondary Infusions Via Gravity Sets:

- Until further notice, **run secondary infusions on a dedicated primary infusion set** instead of via "piggyback" where clinically feasible.
- If infusing via a secondary ("piggyback") line, **always close the upper clamp on the primary infusion set**. Note complete frequent status checks on the infusion to prevent infusion lines from running dry with clamped primary lines.
- Gravity administration sets that do not have clamps above the upper y-site should not be used for secondary infusions.** Attachment 1 includes reference information to provide aid in identifying sets which do not have clamps above the upper Y-site.

2. Secondary Infusions with Infusion Pumps:

- Until further notice, when administering secondary medications via piggyback please **clamp the primary line above the upper Y site using the available clamp** on the pump administration set. This will require unclamping of the primary line after delivery of the secondary medication. The instructions for this differ depending on the type of pump and software as outlined below:

Device	Software Version	Attachment
Infusomat® Space Pump	"U"	2
	"G" and "L"	3
Outlook® ES Pump	All	4
Vista® Basic Pump	All	5

Note: This does not change workflows for rates of infusion greater than 300 mL/hr for Infusomat Pump Administration Sets and 125 mL/hr with Outlook Pump Administration sets where clamping of the primary infusion line is always required per the product IFU to prevent concurrent flow.

- Primary pump administration sets that do not have clamps above the upper y-site should not be used for secondary infusions.** Attachment 1 includes reference information to provide aid in identifying sets which do not have clamps above the upper Y-site.
- Inform all users and post the applicable attachment(s) in care areas where secondary infusions may be administered and where IV administration sets and pump administration sets may be stored until further notice. B. Braun will provide additional notification when users may resume their normal secondary administration process without the need for the modified workflow.
 - Consider priming all IV administration sets and pump administration sets with normal saline or other standard solutions prior to priming with expensive medications to avoid loss of expensive medications or occupational exposure to hazardous medications in the event the set is unable to be primed.
 - If difficulty priming is encountered, obtain a replacement device and file a customer complaint with B. Braun. Adverse reactions or quality problems experienced with the products may be reported to BBMI's Postmarket Surveillance Department by calling 1-833-425-1464.

Actions Required by B. Braun Medical Inc. (BBMI) Customer/User:

- Review this notice in its entirety.



- Ensure that all users in your organization of the above-mentioned products, and other concerned personnel are informed about this voluntary correction and ensure that suitable interim measures are applied as described above. Present information to impacted clinicians at staff meetings, where appropriate.
- Post this notification where the affected products are stored, where secondary infusions may be administered and in break rooms.
- If you are a distributor and have further distributed the product, please forward this notice to your consignees.
- The correction is to be extended to the hospital/healthcare facility level.

2. Return the completed "Urgent Medical Device Correction Acknowledgement Form" to B. Braun Medical Inc. Quality Assurance department by faxing the form to (610) 849-1197 or e-mail to recalls@bbraunusa.com within two (2) weeks of receipt, even if the total inventory in your possession is zero (0).

Adverse Reactions or Quality Problems:

Adverse reactions or quality problems experienced with this product, or questions about this correction may be reported to BBMI's Postmarket Surveillance Department by calling 1-833-425-1464.

In the United States:

Adverse Reactions or quality problems in the United States may also be reported to the FDA's MedWatch Adverse Event Reporting program either online, by regular mail or by fax.

- Complete and submit the report www.fda.gov/medwatch/report.htm
- Regular Mail or Fax: www.fda.gov/MedWatch/getforms.htm or call 1- 800-332-1088 to request a reporting form, then complete and return to the address on the pre-addressed form, or submit by fax to 1-800-FDA-0178

We apologize for the inconvenience that this issue may have caused you and your facility, but we appreciate your understanding of our commitment to assuring our products are safe and effective for both health care professionals and patients. For additional questions regarding this notice, please contact BBMI's Postmarket Surveillance Department at Recalls@bbraunusa.com or by calling 1-844-903-6417.

Sincerely,

Electronically signed by: Cynthia Bauer
Reason: For the reason(s) specified in the document.
Date: Oct 23, 2025 09:11:32 EDT

Cynthia Bauer
Manager II, Customer Complaint Advocacy & Recalls
Postmarket Surveillance
B. Braun Medical Inc.

Enclosures:

Urgent Medical Device Correction Acknowledgement Form

ATTACHMENT 1 – AFFECTED PRODUCTS LIST*

*The products below are subject to the correction. Please note that whether or not the set contains an upper clamp is indicated in column 3. Administration sets that do not have clamps above the upper y-site should not be used for secondary infusions.

Catalog Number	Product Name	Upper Clamp	Primary UDI	Unit of Use-DI	Region Distributed
352049-C	IV ADMIN SET W/ 2 ULTRASITE VLV, 104 IN	Yes	04046964717666	04046964717659	Canada
352062	PB602VSL NON-VENT SET 103-IN-DR*	Yes	04046964181160	04046964181153	United States
352394	US1295 ULTRASITE IV ADMIN SET 60DROP	Yes	04046964293375	04046964293368	United States
352447	V1782F UNIV IV SET W/0.2 MICRON FILTER	No	04046964181665	04046964181658	United States
352604	UNIVERSAL IV ADMIN SET	No	04046964181870	04046964181863	Canada; United States
352635	IV ADMIN SET 60 DROP W/2 SAFEDAY	Yes	04046955594948	04046955594931	United States
352636	IV ADMIN SET 15 DR W/2 SAFEDAY	Yes	04046955594597	04046955594580	United States
352639	IV ADMIN SET 15 DR W/2 SAFEDAY 84 IN.	No	04046955594672	04046955594665	United States
352640	IV ADMIN SET 15 DROP W/3 SAFEDAY	No	04046955594696	04046955594689	United States
352643	IV ADMIN 15D W/2 SAFEDAY & 0.2 FILTER	No	04046955594757	04046955594740	United States
352646	IV ADM SET 15 DR W/2 SAFEDAY & STOPCOCK	No	04046955626946	04046955626939	United States
354202	IV ADMIN SET 15 DROP W/2 CARESITE	No	04046964182037	04046964182020	Canada; Guatemala; United States
354203	IV ADM 15DROP W/2 CARESITE SLIDECLMP	Yes	04046964182051	04046964182044	Canada; United States
354204	IV ADMIN SET 60 DROP W/2 CARESITE	Yes	04046964182075	04046964182068	Canada; United States
354205	IV ADMIN SET 15 DROP W/3 CARESITE	No	04046964182099	04046964182082	Canada; United States
354206	IV ADMIN SET 60 DROP W/3 CARESITE	No	04046964182112	04046964182105	Canada; United States
354207	IV ADM 15 DROP W/2 CARESITE FILTER	No	04046964182136	04046964182129	Canada; United States
354208	IV ADM 15D 2 CARESITE, SAFEPORT MANIFLD	Yes	04046964293818	04046964293801	United States
354209	IV ANES 15DROP W/2 CARESITE ULTRAPORT	No	04046964182150	04046964182143	United States
354210	IV ANES 15DROP W/2 CARESITE 2 ULTRAPORT	No	04046964182174	04046964182167	United States
354211	IV ADM 15 DROP W/3 CARESITE N/C DUAL	Yes	04046964182198	04046964182181	United States
354212	OUTLOOK IV SET 15DROP W/2 CARESITE	Yes	04046964293832	04046964293825	United States
354213	OUTLOOK IV SET 15DROP W/3 CARESITE	Yes	04046964182211	04046964182204	United States
354214	OUTLOOK IV SET 15DR W/2 CARESITE FILTER	Yes	04046964182235	04046964182228	Canada; United States
354302	CARESAFE IV ADMIN SET w/AIRSTOP, 104 IN.	Yes	04046955259946	04046955259939	United States
354303	CARESAFE ADMIN SET w/15um FILTER, 104 IN	Yes	04046955052424	04046955052417	United States
354304	CARESAFE IV ADMIN SET w/AIRSTOP, 114 IN.	No	04046955259960	04046955259953	United States
354305	CARESAFE ADMIN SET w/15um FILTER, 114 IN	No	04046955052448	04046955052431	United States
354306	CARESAFE IV ADMIN SET w/AIRSTOP, 110 IN.	No	04046955259984	04046955259977	United States
354307	CARESAFE ADMIN SET w/15um FILTER, 110 IN	No	04046955052462	04046955052455	United States
354308	CARESAFE ANESTH SET W/AIRSTOP, 134 IN.	No	04046955260003	04046955259991	United States
354310	CARESAFE ANESTH SET w/15um FILTER,134IN.	No	04046955052929	04046955052905	United States
362032	UNIV.15DROP PUMP SET W/3 SAFELINE LL	Yes	04046964182556	04046964182549	Canada; United States

Catalog Number	Product Name	Upper Clamp	Primary UDI	Unit of Use-DI	Region Distributed
362033	UNIV. 15 DR PUMP SET, 1.2 FIL, 2 CRSITE	Yes	04046964837821	04046964837814	Canada; United States
362034	60DROP METRISET PUMP SET, 3 SAFELINE	Yes	04046964182570	04046964182563	United States
362043	Y TYPE BLOOD PUMP SET, SAFELINE INJ	No	04046964182631	04046964182624	United States
362050	UNIV. 15 DROP PUMP SET W/ 0.2 FILTER	Yes	04046964182655	04046964182648	Canada; United States
362420	SPACE PUMP SET 15 DROP W/2 SAFEDAY	Yes	04046955592715	04046955592708	United States
362431	UNIV. 15 DROP PUMP SET W/ 3 ULTRASITE LL	No	04046964294013	04046964294006	United States
362432	SPACE PUMP IV SET, 120 IN.	No	04046964182679	04046964182662	Canada; United States
363010	60 DROP METRISET PUMP SET,3 SFLINE ASV	Yes	04046964294174	04046964294167	Canada; United States
363019	Y-TYPE BLOOD PUMP SET, SAFELINE, ASV	No	04046964294198	04046964294181	United States
363030	UNIV. 15 DROP PUMP SET W/3 SAFELINE, ASV	Yes	04046964294259	04046964294242	United States
363230	UNIV. 15 DROP PUMP SET, 3 ULTRASITE, ASV	Yes	04046964294495	04046964294488	United States
363418	SPACE BLOODSET 200UM FILTER CARESITE ASV	No	04046955856435	04046955856428	United States
363419	Y-TYPE BLOOD PUMP SET, CARESITE, ASV	No	04046964294570	04046964294563	Canada; United States
363420	SPACE PUMP IV SET, 2 CARESITES, ASV	Yes	04046964294594	04046964294587	Canada; United States
363421	UNIV. 15 DROP PUMP SET, 2 CARESITE, ASV	Yes	04046964294617	04046964294600	Canada; United States
363422	60DROP METRISET PUMP, 3 CARESITES, ASV	Yes	04046964294631	04046964294624	Canada; United States
363423	UNIV. 60 DROP PUMP SET, 2 CARESITES, ASV	Yes	04046964294655	04046964294648	Canada; United States
363424	UNIV. 15 DR PUMP SET, 1.2 FIL, 2CRSTE ASV	Yes	04046964837906	04046964837890	Canada; United States
363430	UNIV. 15 DROP PUMP SET, 3 CARESITE, ASV	Yes	04046964294679	04046964294662	Canada; United States
363433	UNIV. 15 DROP PUMP SET, 3 CARESITE, ASV	No	04046964837920	04046964837913	Canada; United States
375100-C	IV ADMIN SET W/ 3 ULTRASITE VLV, 112 IN	No	04046964717680	04046964717673	Canada
375258	RATEFLOW ADMIN 15 DROP W/2 CARESITES	No	04046955256976	04046955256969	United States
456520	ANES SET W/DUAL FEM.LL SAFEPORT MANIFOLD	No	04046964296697	04046964296680	United States
456521	ANESTHESIA SET W/ SAFEPORT MANIFOLD	No	04046964296710	04046964296703	United States
456522	ANES SET W/TRIPLE LL SAFEPORT MANIFOLD	No	04046964296734	04046964296727	United States
457504	IV SET W/2 CARESITE & 2 ULTRAPORT ZERO	No	04046964716089	04046964716072	United States
457506	IV SET W/2 CARESITE & 1 ULTRAPORT ZERO	No	04046964716126	04046964716119	United States
457515	IV SET 3 CARESITE ULTRA ZERO 2/SWAB SAFE	No	04046964716188	04046964716171	United States
470005	IV SET 15DR 3 CARESITE, 124 IN.	No	04046964297199	04046964297182	United States
470179	IV ADMIN SET W/ B/C VALVE,103 IN.	Yes	04046964773181	04046964773174	United States
470183	GLOBAL EXT SET W/CARESITE Y	Yes	04046964719400	04046964715549	Germany; Singapore; United States
473309	AS152V2DLL	No	04046964299636	04046964299629	United States
474912	ET08US EXTENSION SET W/ULTRA Y SITE	Yes	04046964300615	04046964300608	United States
477002	IV ADMIN SET, 15DR, CARESITE LAD, 130 IN	No	04046955202713	04046955202706	United States
477003	IV ADMIN SET, 105 IN.	No	04046955202836	04046955202829	United States
477006	IV ADMIN SET 15D, 4 CARESTIE LAD, 145 IN	Yes	04046955202751	04046955202744	United States
477008	IV ADMIN SET, 121 IN.	No	04046955202850	04046955202843	United States

Catalog Number	Product Name	Upper Clamp	Primary UDI	Unit of Use-DI	Region Distributed
477012	IV ADMIN SET 15D, CARESTIE LAD, 90 IN	Yes	04046955202874	04046955202867	United States
477017	IV ADMIN SET 85 IN.	No	04046955203109	04046955203093	United States
477019	IV ADMIN SET, 2 CARESTIE LAD, 15D, 117IN	Yes	04046955203147	04046955203130	United States
477021	IV ADMIN SET, 15D, 5CARESIE LAD, 136 IN	No	04046955269983	04046955269976	United States
477022	IV ADMIN SET 114 IN.	No	04046955270002	04046955269990	United States
477023	IV ADMIN SET 5 CARESITE LADS, 126 IN.	No	04046955270026	04046955270019	United States
477053	IV SET 15DR 3 SAFEDAY, 106 IN.	No	04046955672431	04046955672424	United States
480231	IV ADMIN SET W/15 DROPS/ML, 128 IN.	No	04046955083534	04046955083428	United States
480233	Y-TYPE BLOOD SET W/ 170µM FILTER	No	04046955083572	04046955083565	United States
480254	SPACE PUMP IV SET 15D, 110 IN.	No	04046955324798	04046955343904	United States
480255	SPACE PUMP IV SET, 127 IN.	Yes	04046955324804	04046955585403	United States
480256	SPACE PUMP SET 15D, CARESITE LAD 122 IN.	Yes	04046955324811	04046955585861	Canada
480258	15 DROP PUMP SET 1 SAFEDAY 0.2 FILTER	Yes	04046955324835	04046955585885	United States
480263	SPACE PUMP IV SET 15D, 90 IN.	No	04046955693931	04046955693924	United States
490025	PB1500 15 DROP IV SET W/2 STOPCOCKS	No	04046964301056	04046964301049	United States
490062	NF1503 15 DR UNIV SPIKE, INJ SITE, BCV	No	04046964301445	04046964301438	United States
490067	IV SET 15DR 5CARESIE 1 ULRAPORT, 141 IN	No	04046964186042	04046964186035	United States
490070	NF1500UP, IV ADMINISTRATION SET	No	04046964301582	04046964301575	United States
490088	NF1505, IV ADMIN SET W/2 H.F. STOPCOCKS	No	04046964301889	04046964301872	United States
490100	UNIV. 15 DROP PUMP SET W 3 CARESITE LL	Yes	04046964186080	04046964186073	Canada; United States
490101	60 DROP METRISET PUMP SET, 3 CARESITES	Yes	04046964186103	04046964186097	Canada; United States
490102	UNIV. 15 DROP PUMP SET W/2 CARESITE CV	Yes	04046964186127	04046964186110	Canada; Germany; United States
490103	SPACE PUMP IV SET W/2 CARESITE, CKVLV	Yes	04046964302107	04046964302091	Canada; United States
490104	UNIV 60 DROP PUMP SET, W/2 CARESITE INJ.	Yes	04046964186141	04046964186134	Canada; United States
490105	Y-TYPE BLOOD PUMP SET, CARESITE INJ	No	04046964186165	04046964186158	Canada; United States
490107	IV ADMIN SET W/2 SAFELINE & HF STOPCOCK	No	04046964302145	04046964302138	United States
490159	BURETTE SET W/ 2 CARESITE, 134 IN.	Yes	04046964303128	04046964303111	United States
490189	IV ADMIN SET W/15 DROPS/ML, 108 IN.	No	04046964303586	04046964303579	United States
490193	IV ADMIN SET 15DR HF 4-WAY STPCK 135 IN.	Yes	04046964303661	04046964303654	United States
490196	Y-TYPE BLOOD SET 124 IN.	No	04046964303722	04046964303715	United States
490202	EXT SET W/ B/C VLV, ULTRAPORT, 57 IN.	No	04046964303845	04046964303838	United States
490208	IV ADMIN SET W/15 DR / mL W/ EXT SET	Yes	04046964303968	04046964303951	United States
490225	IV ADMIN SET, 129 IN, 3 CARESITE LADS	Yes	04046964304248	04046964304231	United States
490226	IV ADMIN SET 137 IN.	No	04046964304262	04046964304255	United States
490234	EXT SET ULTRAPORT STOPCOCK, 37IN.	No	04046964716362	04046964716355	United States
490238	EXT SET W/ 6 B/C VALVES, 17 IN	No	04046964716409	04046964716393	United States
490239	IV ADMIN SET 103 IN.	No	04046964716423	04046964716416	United States

Catalog Number	Product Name	Upper Clamp	Primary UDI	Unit of Use-DI	Region Distributed
490241	IV ADMIN SET HF STOPCOCK,107 IN.	No	04046964716447	04046964716430	United States
490242	IV ADMIN SET 60 DROP/ML, 105 IN.	No	04046964716461	04046964716454	United States
490245	IV ADMIN SET W/15 DROPS/ML, 133 IN.	No	04046964716508	04046964716492	United States
490252	BURETTE SET, 4 CARESITE LADS, STPCK 125"	Yes	04046955545636	04046955545629	United States
490274	IV ADMIN SET 3 CARESITE LADS, 131 IN.	No	04046964716621	04046964716614	United States
490283	IV ADMIN SET W/15 DROPS/ML 126 IN	No	04046964716669	04046964716652	United States
490291	EXTENSION SET W/ ULTRAPORT, 48 IN.	No	04046964716768	04046964716751	United States
490298	IV ADMIN SET W/15 DROPS/ML, 116 IN.	No	04046964716829	04046964716812	United States
490301	IV AD.SET W/15 DR 128 IN, EXT SET 9 IN	No	04046964716867	04046964716850	United States
490303	IV ADMIN SET,3 CARESITE LADS, 137 IN	No	04046964526381	04046964526374	United States
490309	IV ADMIN SET W/ 2 CARESITE, 140 IN	No	04046964527517	04046964527500	United States
490311	IV ADMN SET W/15 DROP,2 CARESITE,142 IN.	Yes	04046964539688	04046964525629	United States
490312	IV ADMN SET W/15 DROP, ULTRAPORT,134 IN.	No	04046964525643	04046964525636	United States
490315	IV ADMIN SET, 15 DR/ML,2 CARESITE,146 IN	No	04046964527173	04046964525612	United States
490317	IV AD SET CARESITE 127 IN. EXT SET 9 IN.	No	04046964525742	04046964525735	United States
490326	IV AD SET, 2 CARESITE, SAFEPORT, 134 IN.	No	04046964541063	04046964541056	United States
490329	IV ADMIN SET, 81 IN W/ 32 IN EXT SET	No	04046964541100	04046964541094	United States
490331	IV ADMN SET,2 CARESITE, SAFEPORT,129 IN.	No	04046964541148	04046964541131	United States
490340	IV ADMIN SET, QUAD MANIFOLD,120 IN.	No	04046964590481	04046964590474	United States
490343	IV ADMIN SET,2 CARESITE, TRIP MANIFOLD	No	04046964593659	04046964593642	United States
490349	IV ADMIN SET W/ 60 DROPS/ML,141 in.	No	04046964602405	04046964602399	United States
490353	IV ADMIN SET W/15 DROPS/ML 115 in.	No	04046964602368	04046964602351	United States
490354	IV ADMIN SET W/15 DROPS/ML 115 in.	No	04046964602382	04046964602375	United States
490355	Y-TYPE BLOOD SET W/ 170 M FILTER 120 IN	Yes	04046964616013	04046964616006	United States
490357	IV ADMIN SET W/ 10 DROPS/ML, 134 in.	No	04046964644399	04046964644382	United States
490358	EXTENSION SET W/ B/C VALVE 53 in.	No	04046964642746	04046964642739	United States
490360	IV ADMIN SET,15 DROP,4 CARESITE, 119 IN.	No	04046964642609	04046964642487	United States
490361	OUTLOOK PUMP SET,3 CARESITE LADS,115 IN.	Yes	04046964642623	04046964642616	United States
490367	IV ADMIN SET W/60 DR/ML, CARESITE,147 IN.	No	04046964669903	04046964669798	United States
490376	IV ADMIN SET 15 DROPS/ML, CARESITE,98 IN	No	04046964718182	04046964718175	United States
490385	IV ADMIN SET, 10 DROPS/ML, 144 IN	No	04046964730146	04046964730139	United States
490389	IV ADMIN SET 15 DR/ML, 2 B/C VLV, 97IN.	No	04046964730207	04046964730191	United States
490392	IV ADMIN SET W/60 DROPS/ML 106 IN.	No	04046964737367	04046964737350	United States
490400	IV ADMIN SET W/CARESITE EXT SET	No	04046964761393	04046964761386	United States
490401	IV ADMIN SET, 3 CARESITES, 120 IN.	No	04046964761416	04046964761409	United States
490407	IV ADMIN SET W/CARESITE EXT SET	Yes	04046964773242	04046964773235	United States
490408	IV ADMIN SET , 2 CARESITES, 125 IN.	No	04046964773266	04046964773259	United States

Catalog Number	Product Name	Upper Clamp	Primary UDI	Unit of Use-DI	Region Distributed
490412	130 IN.IV ADMIN SET W/ 9 IN. EXT SET	No	04046964797989	04046964797972	United States
490426	IV ADMIN SET W/10 DROPS, 138 IN.	No	04046964823503	04046964823497	United States
490427	IV ADMIN SET W/ 15 DROPS/ML 124 IN.	No	04046964823527	04046964823510	United States
490428	IV ADMIN SET W/ 60 DROPS/ML 120 IN.	No	04046964823541	04046964823534	United States
490432	IV ADMIN SET W/ 15 DROPS/ML, 120 IN.	No	04046964832826	04046964832819	United States
490437	IV ADMIN SET W/ 15 DROPS/ML, 118 IN.	No	04046964833427	04046964833410	United States
490438	Y-TYPE BLOOD SET, HAND PUMP, 146 IN.	No	04046964833441	04046964833434	United States
490444	IV ADMIN SET, 15D 2 CARESITE LAD, 126 IN	No	04046964848773	04046964848766	United States
490448	EXTENSION SET, 18 IN	Yes	04046964848858	04046964848841	United States
490451	IV ADMIN SET 124 IN.	No	04046964848919	04046964848902	United States
490452	IV ADMIN SET, 15 DROPS/ML 125 IN	No	04046964848933	04046964848926	United States
490465	EXTENSION SET, 52 IN.	No	04046964950735	04046964950728	United States
490467	IV AD. SET,15 DROPS/ML, W/EXT SET	No	04046964958779	04046964958762	United States
490468	IV ADMIN SET 117 IN.	No	04046964958793	04046964958786	United States
490469	IV ADMIN SET W/15 DROPS/ML, 132 IN	No	04046964958816	04046964958809	United States
490470	IV ADMIN SET, 60D 2 CARESITE VLV, 132 IN	No	04046964958830	04046964958823	United States
490472	IV ADMIN SET, 15D, 144 IN.	No	04046964958878	04046964958861	United States
490476	IV ADMIN SET, 2-HIGH FLOW STPKCS, 124 IN	No	04046964958953	04046964958946	United States
490477	IV ADMIN SET W/15 DROPS/ML, 134 IN.	No	04046964991271	04046964991264	United States
490484	IV ADMIN SET W/15 DROPS, 162 IN.	No	04046964991219	04046964991202	United States
490486	EXTENSION SET 49 IN.	No	04046964991257	04046964991240	United States
490487	IV ADMIN SET W/15 DROPS/ML, 125 IN.	No	04046955086856	04046955086849	United States
490488	IV ADMIN SET W/15 DROPS/ML, 126 IN.	No	04046955086870	04046955086863	United States
490489	IV ADMIN SET 15DR, 92 IN W/EXT SET 9 IN.	No	04046955086894	04046955086887	United States
490491	IV ADMINISTRATION SET, 138 IN.	No	04046955086931	04046955086924	United States
490492	IV ADMIN SET W/15 DROPS/ML, 132 IN.	No	04046955086757	04046955086740	United States
490503	IV ADMIN SET W/15 DROPS/ML, 139 IN.	No	04046955135714	04046955135707	United States
490504	IV ADMIN SET W/15 DROPS/ML, 108 IN.	No	04046955135738	04046955135721	United States
490507	IV ADMINISTRATION SET WITH 15 DROPS/ML	No	04046955135790	04046955135783	United States
490508	IV ADMIN SET W/ 15DROPS/ML	No	04046955135813	04046955135806	United States
490511	IV ADMIN SET W/ 15 DROPS/ML, 99 IN	No	04046955135875	04046955135868	United States
490513	IV ADMIN SET, 15DR, CARESITE LAD, 121 IN	No	04046955137237	04046955137220	United States
490519	IV ADMIN SET 15 DROP, 132 IN	No	04046955185238	04046955186549	United States
490521	IV ADMIN SET, 15D, 120 IN.	No	04046955185252	04046955186563	United States
490522	IV ADMIN SET, 15DR, CARESITE LAD, 116 IN	No	04046955203123	04046955203116	United States
490527	IV ADMIN SET 15D 2 CARESITE LADS, 112IN.	No	04046955203222	04046955203215	United States
490528	IV ADMIN SET 15 DROP W/ CARESITE EXT SET	No	04046955203284	04046955203277	United States

Catalog Number	Product Name	Upper Clamp	Primary UDI	Unit of Use-DI	Region Distributed
490529	IV ADMIN SET 15 DRP CARESITE SAFEPORT	No	04046955203307	04046955203291	United States
490530	Y-TYPE BLOOD SET, 3 CARESITE LADS, 146 IN	No	04046955203321	04046955203314	United States
490535	IV ADMIN SET, 15D, 139 IN.	No	04046955203406	04046955203390	United States
490542	IV ADMIN SET 15D, 2 CARESITES, 153 IN.	No	04046955270187	04046955270170	United States
490543	IV ADMIN SET 15D, CARESITE LAD, 122IN.	Yes	04046955270200	04046955270194	United States
490546	IV ADMIN SET 15D 3 CARESITE LADS, 150IN.	No	04046955270262	04046955270255	United States
490547	IV ADMIN SET, 125 IN.	No	04046955270286	04046955270279	United States
490549	IV ADMIN SET, 150 IN.	No	04046955270323	04046955270316	United States
490550	111 IN. IV ADMIN SET W/ 12 IN. Y-EXT SET	No	04046955270347	04046955270330	United States
490551	AESTHESIA IV SET, 127 in.	No	04046955275304	04046955275298	United States
490552	IV ADMIN SET 15D, 5 CARESITE INJ, 150in.	No	04046955275342	04046955275335	United States
490553	IV ADMIN SET, 120IN.	No	04046955275366	04046955275359	United States
490554	IV ADMIN SETW/ B/C VLV 3 CARESITE 134IN.	Yes	04046955275380	04046955275373	United States
490555	IV ADMIN SET 15DR w/3 B/C VLV, 124 IN.	No	04046955275403	04046955275397	United States
490558	IV ADMIN SET 15D w/B/C VLV, 139 IN.	No	04046955275465	04046955275458	United States
490559	IV ADMINISTRATION SET, 117 IN.	No	04046955275489	04046955275472	United States
490560	IV ADMINISTRATION SET, 108 IN.	No	04046955275502	04046955275496	United States
490561	IV ADMINISTRATION SET, 134 IN.	No	04046955275526	04046955275519	United States
490562	IV ADMIN SET 15D w/B/C VLV, 140 IN.	No	04046955275540	04046955275533	United States
490563	IV ADMIN SET, 111 IN.	No	04046955275564	04046955275557	United States
490564	OUTLOOK PUMP SET 3 CARESITE LADS 133 IN.	Yes	04046955275588	04046955275571	United States
490565	IV ADMIN SET, 91 IN.	No	04046955275601	04046955275595	United States
490566	IV ADMIN SET, 110 IN.	No	04046955275625	04046955275618	United States
490568	IV ADMIN SET 115, IN. W. 9 IN EXT SET	No	04046955275663	04046955275656	United States
490569	IV ADMIN SET 15D 2 CARESITE LADS, 122 IN	No	04046955275687	04046955275670	United States
490573	IV ADMIN SET 15D STOPCOCK, DUAL MANIFOLD	No	04046955287451	04046955287444	United States
490574	IV ADMIN SET 15D w/B/C VLV, 110 IN.	No	04046955287475	04046955287468	United States
490575	IV ADMIN 15D 5 CARESITE LADS STPCK 120in	No	04046955287499	04046955287482	United States
490577	EXTENSION SET, 12 CARESITE LADS	Yes	04046955287536	04046955287529	United States
490578	IV ADMIN SET 15D 5 CARESITE LADS, 140 IN	No	04046955287550	04046955287543	United States
490581	IV ADMIN SET W/15D, 4 CARESITE LADS	No	04046955287611	04046955287604	United States
490582	IV ADMIN 15D CARESITE LADS, STPCK 134 IN	No	04046955287635	04046955287628	United States
490583	IV ADMIN SET, 15D, 108 IN.	Yes	04046955287659	04046955287642	United States
490588	IV ADMIN SET, 10 DROP, 130 IN.	No	04046955287758	04046955287741	United States
490589	IV ADMIN SET, 132 IN.	No	04046955287772	04046955287765	United States
490600	IV ADMIN SET, 15DR, 105 IN.	No	04046955352777	04046955352760	United States
490605	IV ADMIN SET, 125 IN. w/ EXT SET 9 IN.	No	04046955352876	04046955352869	United States

Catalog Number	Product Name	Upper Clamp	Primary UDI	Unit of Use-DI	Region Distributed
490607	IV ADMIN SET,3 CARESITE,3 HF SPCK, 129IN	No	04046955352913	04046955352906	United States
490608	IV SET 15DR 3CARESITE 2STOPCOCK, 126IN	No	04046955352937	04046955352920	United States
490612	15DR 2SAFELINE 1CARESITE 1STOPCOCK 131IN	No	04046955664610	04046955664603	United States
490615	IV ADMIN SET 15D 2 CARESITE LADS, 138 IN	No	04046955664672	04046955664665	United States
490616	IV SET 15 DR 3CARESITE 2 ULTRAPORT	Yes	04046955664696	04046955664689	United States
490618	IV SET 15DR 2SAFEDAY 2SAFEPORT, 99 IN.	No	04046955664733	04046955664726	United States
490619	IV SET 60DR 2SAFEDAY 2SAFEPORT, 99 IN.	No	04046955664757	04046955664740	United States
490629	EXTENSION SET, 25 IN.	No	04046955665051	04046955665044	United States
490644	IV ADMIN SET 15 DRPS 4W STOPCOCK W/ C&S	No	04046955765003	04046955764990	United States
490653	IV ADMIN 15D, 2 B/C VLV, 3 CARESITE LADS	No	04046955836338	04046955836345	United States
490654	IV ADMIN SET, 2 CARSITE LADS	No	04046955836369	04046955836352	United States
490655	UNV 15 DR PMP ST W/3 CARESITE LL 113.2	Yes	04046955836383	04046955836376	United States
490657	IV ADMIN SET 15 DRPS ULTRAP. ZERO 90.6IN	No	04046955836420	04046955836413	United States
490658	IV ADMIN SET 15 DRPS STOPCOCK MALE 117IN	No	04046955836697	04046955836680	United States
NF1251	SAFELINE ADDITIV ALTVTD 3 INJ S	No	04046964187100	04046964187094	Canada; United States
NF1270A	NF1270A SAFELINE ANESTHESIA IV SET	No	04046964187124	04046964187117	United States
NF1290	SAFELINE PED. ADDITIVE	No	04046964187148	04046964187131	United States
US9923F	VISTA BASIC PUMP IV SET W/ 15 DROP/ML	No	04046964838156	04046964838149	United States
V1443	ADULT ADD SET CHK VLV 2 INJ S	No	04046964188350	04046964188343	Canada; United States
V1447	ADDITIVE SET W/ 3 INJ SITES -DR*	No	04046964188411	04046964188404	Canada
V1448	ADD PEDIATRIC 3 INJ LUER LOCK	No	04046964188435	04046964188428	Canada; United States
V1484	ADULT ADD 2 INJ CKV LUER LOCK	No	04046964188473	04046964188466	United States
V5405	EXT/2INJ/FLBL/FML LRLK/CV 34	Yes	04046964188978	04046964188961	United States

ATTACHMENT 2

SECONDARY ADMINISTRATION PROCESS INFUSOMAT SPACE® SOFTWARE U

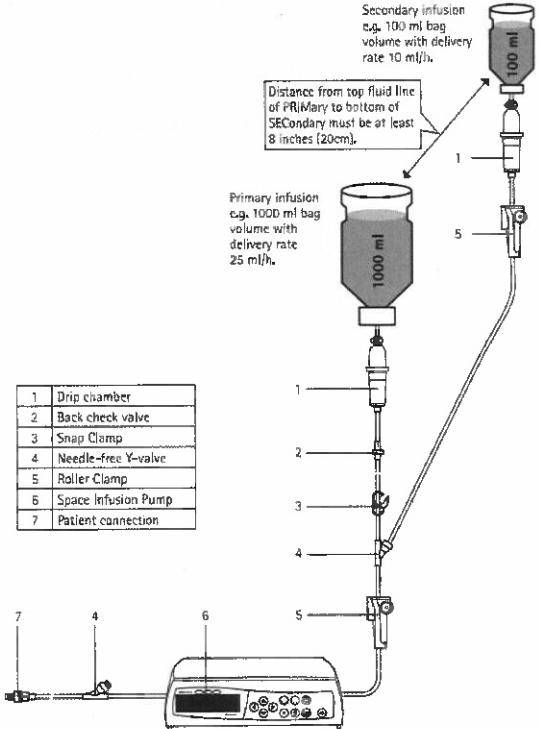
PURPOSE:

To provide workflow recommendations for clamping of the primary infusion line when administering secondary medications.

APPLICABILITY:

These recommendations apply only for pump administration sets used within Infusomat® Space Pumps with Software "U" versions. The software version of the Infusomat Space Pump can be verified by accessing the pump "Options" -> "Status" menu. Menus are accessed from the Home screen using the  arrow keys. Press  key while pump is running to access the Home screen. Software "U" versions begin with the prefix "X86U" (Example: X86U030017).

Note: Primary infusions with no secondaries connected can continue to be infused without impact. Refer to product IFU.

Step	Instruction
1	Close roller clamp. Fully spike secondary container. Hang fluid container. Gently squeeze and release drip chamber of secondary line until approximately half full. <i>NOTE: Vent cap should remain in the closed position unless infusing from a solid stopper glass bottle.</i> <i>NOTE: SECondary infusion should not be used for critical infusions.</i>
2	Prime secondary set by slowly opening roller clamp and allow fluid to fill tubing. Close roller clamp. Ensure air is expelled. Repeat prime if necessary.
3	Lower primary container such that the top fluid line of the primary is at least 8 in. (20 cm) below the bottom of the secondary IV container. 
4	Swab upper access port on primary set and attach secondary set to the primary set.

B | BRAUN

Step	Instruction	
5	Close the clamp on the primary administration set above the upper Y-Site.	
6	Follow the product IFU to enter the Secondary menu. Arrow down to the setting "Back to PRIM", press the left arrow key, and change the setting to "Manual". This will cause the pump to alarm when the secondary administration is completed. Note: The method of switching back to PRIMary remains for all SECondary infusions until changed. Program Secondary infusion per IFU. Note: If you do not see this option, decrease the Primary VTBI to less than 1mL before programming Secondary.	
7	Allow pump to infuse SECondary. When the VTBI reaches 0 the pump will alarm; to notify the clinician that the SECondary is complete. Ensure volume is complete. If not, program additional volume. Verify all infusion data is as intended prior to restarting. If volume is complete, unclamp the clamp on the PRIMary line and close the roller clamp on the SECondary line.	
8	Arrow up on SECondary home page -> use the left arrow to Resume PRIMary and press Start to infuse the PRIMary line.	

ATTACHMENT 3

SECONDARY ADMINISTRATION PROCESS INFUSOMAT SPACE® SOFTWARE G & L

PURPOSE:

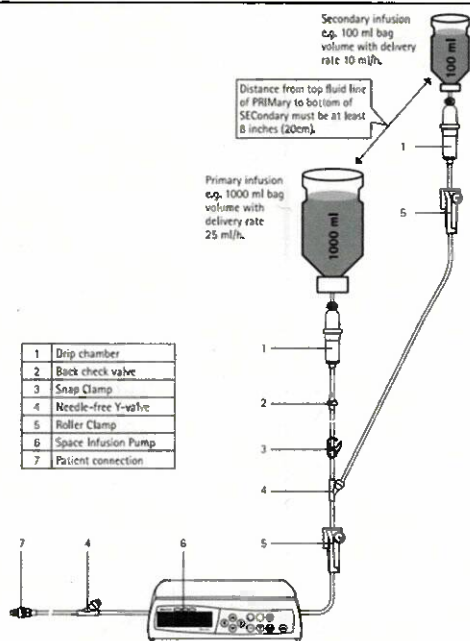
To provide workflow recommendations for clamping of the primary infusion line when administering secondary medications.

APPLICABILITY:

These recommendations apply only for pump administration sets used within Infusomat® Space Pumps with Software "G" or "L" versions. The software version of the Infusomat Space Pump can be verified by accessing the pump "Options" -> "Status" menu. If your device does not allow for this workflow please contact B. Braun Medical at medinfo.us@bbraunusa.com or call 1800-854-6851 as settings within the configuration file may need to be changed. Menus are accessed from the Home screen using the  arrow keys. Press  key while pump is running to access the Home screen. Software "G" versions begin with the prefix "X86G" (Example: X86G030102). Software "L" versions begin with the prefix "X86L" (Example: X86L032101). Note that Software L is utilized only in Canada.

Note: Primary infusions with no secondaries connected can continue to be infused without impact. Refer to product Instructions for Use (IFU).

Step	Instruction
1	Close roller clamp. Fully spike secondary container. Hang fluid container. Gently squeeze and release drip chamber of secondary line until approximately half full. <i>NOTE: Vent cap should remain in the closed position unless infusing from a solid stopper glass bottle.</i> <i>NOTE: SECondary infusion should not be used for critical infusions.</i>
2	Prime secondary set by slowly opening roller clamp and allow fluid to fill tubing. Close roller clamp. Ensure air is expelled. Repeat prime if necessary.
3	Lower primary container such that the top fluid line of the primary is at least 8 in. (20 cm) below the bottom of the secondary IV container.



4 Swab upper access port on primary set and attach secondary set to the primary set.

5 Close the clamp on the primary administration set above the upper Y-Site.



6 Go to program the Piggyback through Special Functions ->Piggyback->enter drug information as need be. When prompted, to select automatic change yes or no, choose "no" for a manual change.
Note: If you do not see this option, decrease the Primary VTBI to less than 1mL before programming Secondary.

7 Program the Piggyback infusion with appropriate clinical data specific to the needs of the patient and treatment being delivered (e.g. infusion rate, Volume To Be Infused, etc.), noting the VTBI to be as accurate as possible.

8 Allow pump to infuse Piggyback. When the VTBI reaches 0 the pump will alarm, to notify the clinician that the Piggyback is complete. Ensure volume is complete. If not, program additional volume. Verify all infusion data is as intended prior to restarting. If volume is complete, unclamp the clamp on the PRIMARY line and close the roller clamp on the Piggyback line.



9 Arrow up on Piggyback home page -> use the left arrow to Resume PRIMARY and press Start to infuse the PRIMARY line.

ATTACHMENT 4

SECONDARY ADMINISTRATION PROCESS OUTLOOK ES®

PURPOSE:

To provide workflow recommendations for clamping of the primary infusion line when administering secondary medications.

APPLICABILITY:

These recommendations apply only for pump administration sets used within Outlook ES® pumps.

Note: Primary infusions with no secondaries connected can continue to be infused without impact. Refer to product Instructions for Use (IFU).

Step	Instruction	Picture
1	Close roller clamp. Fully spike secondary container. Hang fluid container. Gently squeeze and release drip chamber of secondary line until approximately half full. <i>NOTE: Vent cap should remain in the closed position unless infusing from a solid stopper glass bottle.</i> <i>NOTE: SECondary infusion should not be used for critical infusions.</i>	
2	Prime secondary set by slowly opening roller clamp and allow fluid to fill tubing. Close roller clamp. Ensure air is expelled. Repeat prime if necessary.	
3	Lower primary container such that the top fluid line of the primary is at least 8 in. (20 cm) below the bottom of the secondary IV container.	Secondary infusion e.g. 100 ml bag volume with delivery rate 10 ml/h. Distance from top fluid line of PRIMARY to bottom of SECONDARY must be at least 8 inches (20cm). Primary infusion e.g. 1000 ml bag volume with delivery rate 25 ml/h. 1 Drip chamber 2 Back check valve 3 Snap Clamp 4 Needle-free Y-valve 5 Roller Clamp 6 Infusion Pump 7 Patient connection Outlook Pump
4	Swab upper access port on primary set and attach secondary set to the primary set.	

B | BRAUN

5	Close the clamp on the primary administration set above the upper Y-Site.	
6	Ensure the pump is programmed to "Piggyback Callback" ON. In this mode, the pump will alarm at the end of the SECondary Infusion and the PRIMary will not automatically begin to infuse. To move the pump to "Piggyback call back ON" from the Home Screen- Press the , -> then use  to arrow to "Piggyback Callback" and press Enter. Cursor over to ON with the screen. Note: This setting will return to OFF when the pump is turned off.	Piggyback Callback When Piggyback Callback is ON, KVO will alarm as soon as the Piggyback is done. Set Callback: OFF  Cursor  & Press ENTER
7	Program the SECondary infusion with appropriate clinical data specific to the needs of the patient and the treatment being delivered (e.g. infusion rate, Volume To Be Infused, etc.), noting the VTBI to be as accurate as possible.	
8	Allow pump to infuse SECondary. When the VTBI reaches 0 the pump will alarm; to notify the clinician that the piggyback is complete.	
9	Ensure volume is complete. If not, program additional volume. Note: Verify all infusion data is as intended prior to restarting. If volume is complete, unclamp the clamp on the primary line and close the roller clamp on the piggyback/secondary line.	

ATTACHMENT 5

SECONDARY ADMINISTRATION PROCESS VISTA® BASIC

PURPOSE:

To provide workflow recommendations for clamping of the primary infusion line when administering secondary medications.

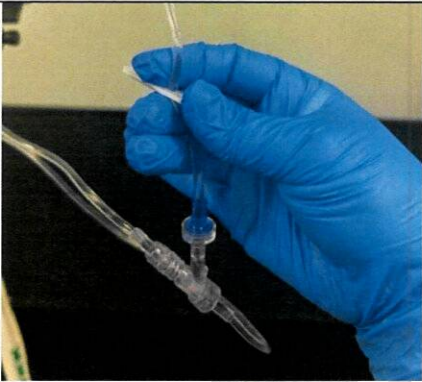
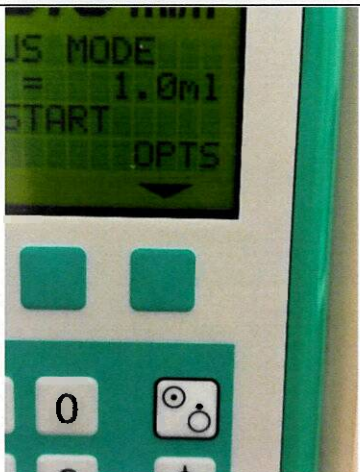
APPLICABILITY:

These recommendations apply only for pump administration sets used within VISTA® Basic pumps.

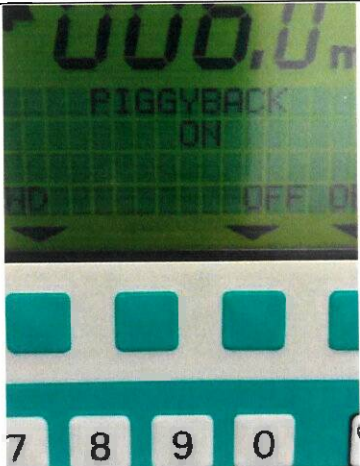
Note: Primary infusions with no secondaries connected can continue to be infused without impact. Refer to product Instructions for Use (IFU).

Step	Instruction
1	Close roller clamp. Fully spike secondary container. Hang fluid container. Gently squeeze and release drip chamber of secondary line until approximately half full. <i>NOTE: Vent cap should remain in the closed position unless infusing from a solid stopper glass bottle.</i> <i>NOTE: SECondary infusion should not be used for critical infusions.</i>
2	Prime secondary set by slowly opening roller clamp and allow fluid to fill tubing. Close roller clamp. Ensure air is expelled. Repeat prime if necessary.
3	Lower primary container such that the top fluid line of the primary is at least 8 in. (20 cm) below the bottom of the secondary IV container. Secondary infusion e.g. 100 ml bag volume with delivery rate 10 ml/h. Distance from top fluid line of PRIMARY to bottom of SECondary must be at least 8 inches (20cm). Primary infusion e.g. 1000 ml bag volume with delivery rate 25 ml/h. 1 Drip chamber 2 Back check valve 3 Snap Clamp 4 Needle-free Y-valve 5 Roller Clamp 6 Infusion Pump 7 Patient connection Vista Pump
4	Swab upper access port on primary set and attach secondary set to the primary set.

B | BRAUN

5	Close the clamp on the primary administration set above the upper Y-Site.	
6	The pump automatically switches from the Piggyback to the Continuous Mode infusion when the Piggyback program is complete. If you program the infusion as a Piggyback, ensure that your Primary volume is programmed to less than 1mL. This will ensure the line does not run dry and the pump will alarm after the secondary infusion is complete. Program the Piggyback infusion with appropriate clinical data, noting the VTBI to be as accurate as possible.	
7	Allow pump to infuse Piggyback. When the VTBI reaches 0 the pump will alarm, to notify the clinician that the Piggyback is complete.	
8	Press END.	
9	Press OPTS.	

B | BRAUN

10	Press OFF to turn Piggyback off, open clamp on Primary tubing, and close the roller clamp on the Piggyback line.	
11	Press END to return to Continuous Mode.	



URGENT MEDICAL DEVICE CORRECTION ACKNOWLEDGMENT FORM

IV ADMINISTRATION SETS & PUMP ADMINISTRATION SETS (INFUSOMAT® SPACE; OUTLOOK®; & VISTA BASIC) Occlusion and Backflow

Please complete and return this form by one of the below options:

Fax: (610) 849-1197 or Email recalls@bbraunusa.com

Acknowledgement of Correction

Completing and returning this form indicates you have read and understand the Urgent Medical Device Correction Notification. If you have further distributed these impacted products, please forward the Urgent Medical Device Correction Notification Letter to your affected consignees.

☐ ***DISTRIBUTORS ONLY*** – checking this box indicates you or your Corporate Office have forwarded this notice to all affected customers.

PLEASE FULLY COMPLETE:

Customer Information	
Facility Name: The Mount Sinai Hospital Of Queens	Printed Name:
Address: 30-14 Crescent St, LONG ISLAND CITY, NY 11102	Signature of Person Completing Form:
Account Number: 0020073586	Telephone Number:
Fax Number:	Date:
E-Mail:	

