



Account number: 6213735

CAPA 2025-007

LeMaitre Vascular, Inc.

24 APR 2025

URGENT: MEDICAL DEVICE RECALL**TufTex® Over-the-Wire Embolectomy Catheter****Pruitt® Irrigation Occlusion Catheter****Pruitt® Occlusion Catheter****(1) Attention to Customer:**

Attn: Risk Management Department
MOUNT SINAI SOUTH NASSAU
2445 OCEANSIDE ROAD
OCEANSIDE, NY 11572
UNITED STATES OF AMERICA

Dear Device Customer/Distributor,

(2) Purpose of this letter

The purpose of this letter is to advise you that LeMaitre Vascular, Inc. is voluntarily recalling the following devices.

Product Line	Catalog Numbers	Lot Numbers
TufTex® Over-the-Wire Embolectomy Catheter	1651-34, 1651-38, 1651-44, 1651-48, 1651-68, 1651-78, 1651-84, 1651-88, E1651-34, E1651-38, E1651-44, E1651-48, E1651-64, E1651-68, E1651-78, E1651-84, E1651-88	OTW4001, OTW4176, OTW4177, OTW4181, OTW4183, OTW4184, OTW4189, OTW4190, OTW4191, OTW4192, OTW4193, OTW4194, OTW4195, OTW4196, OTW4211, OTW4212, OTW4231, OTW4232, OTW4233, OTW4266, OTW4267, OTW4268, OTW4269, OTW4270, OTW4271, OTW4272, OTW4273, OTW4274, OTW4275, OTW4276, OTW4277, OTW4283, OTW4284, OTW4285, OTW4286, OTW4287, OTW4288, OTW4292, OTW4293, OTW4300, OTW4301, OTW4302, OTW4305, OTW4311, OTW4317, OTW4319, OTW4320, OTW4328, OTW4333, OTW4334, OTW4336, OTW4346, OTW4347, OTW4348, OTW4349, OTW4360, OTW4362, OTW4366, OTW4368, OTW4377, OTW4384, OTW4387, OTW4388, OTW4389, OTW4391, OTW4392, OTW4395, OTW4396, OTW4397, OTW4398, OTW4399, OTW4400, OTW4405, OTW4406, OTW4407, OTW4417, OTW4429, OTW4430, OTW4431, OTW4432, OTW4433, OTW4436, OTW4437, OTW4444, OTW4445, OTW4447, OTW4452, OTW4453, OTW4454, OTW4456, OTW4457, OTW4458, OTW4459A, OTW4461, OTW4462, OTW4469, OTW4470, OTW4471, OTW4472, OTW4473, OTW4474, OTW4475, OTW4477, OTW4480, OTW4481, OTW4482, OTW4483, OTW4484, OTW4486, OTW4487, OTW4488, OTW4491, OTW4492, OTW4493, OTW4502, OTW4507, OTW4508, OTW4512, OTW4513, OTW4515, OTW4517, OTW4518, OTW4522A, OTW4523, OTW4524, OTW4525, OTW4526, OTW4528, OTW4529A, OTW4533, OTW4534, OTW4535, OTW4536, OTW4544, OTW4548, OTW4555, OTW4558, OTW4560, OTW4561, OTW4576, OTW4577, OTW4578, OTW4579, OTW4580, OTW4581, OTW4584, OTW4588, OTW4602, OTW4614, OTW4619, OTW4620, OTW4621, OTW4622, OTW4623, OTW4624, OTW4644, OTW4646, OTW4648, OTW4661, OTW4664, OTW4674, OTW4686, OTW4687, OTW4688, OTW4696, OTW4699, OTW4702, OTW4707, OTW4709, OTW4716, OTW4721, OTW4723, OTW4724, OTW4735, OTW4775, OTW4794, QOT1001, QOT1002, QOT1003, QOT1004, QOT1005, QOT1006, QOT1007, QOT1008, QOT1009, QOT1010, QOT1012, QOT1013, QOT1014, QOT1015, QOT1016, QOT1017, QOT1018, QOT1019, QOT1020, QOT1022, QOT1023, QOT1024, QOT1025, QOT1031, QOT1035, QOT1040, QOT1048, QOT1051, QOT1052, QOT1053, QOT1061, QOT1062, QOT1087, QOT1088, QOT1089, QOT1095, QOT1097, QOT1099, QOT1101, QOT1105, QOT1108, QOT1109, QOT1110, QOT1116, QOT1117, QOT1118, QOT1120, QOT1122, QOT1123, QOT1128, QOT1133, QOT1134, QOT1135, QOT1139, QOT1145, QOT1148, QOT1150, QOT1151, QOT1152, QOT1153, QOT1154, QOT1161, QOT1162,



Account number: 6213735

CAPA 2025-007

Product Line	Catalog Numbers	Lot Numbers
		QOT1163, QOT1169, QOT1170, QOT1171, QOT1175, QOT1176, QOT1177, QOT1178, QOT1181, QOT1193, QOT1206, QOT1207, QOT1209, QOT1210, QOT1211, QOT1212, QOT1222, QOT1225, QOT1236, QOT1237, QOT1243, QOT1244, QOT1245, QOT1247, QOT1250, QOT1255, QOT1256, QOT1258, QOT1266, QOT1267, QOT1268, QOT1271, QOT1273, QOT1277, QOT1278, QOT1287, QOT1288, QOT1289, QOT1291, QOT1292, QOT1302, QOT1303, QOT1308, QOT1310, QOT1311, QOT1317, QOT1322, QOT1325, QOT1330, QOT1333, QOT1334, QOT1336, QOT1338, QOT1340, QOT1342, QOT1344, QOT1347, QOT1355, QOT1357, QOT1358, QOT1359, QOT1366, QOT1367, QOT1368, QOT1372, QOT1374, QOT1375, QOT1381, QOT1382, QOT1383, QOT1384, QOT1387, QOT1389, QOT1394, QOT1399, QOT1406, QOT1407, QOT1408, QOT1409, QOT1412, QOT1416, QOT1420, QOT1424, QOT1425, QOT1430, QOT1439, QOT1441, QOT1443, QOT1444, QOT1445, QOT1467, QOT1468, QOT1473, QOT1483, QOT1485, QOT1488, QOT1490, XOT00001601, XOT00001614, XOT00001615, XOT00001616, XOT00001617, XOT00001618, XOT00001634, XOT00001635, XOT00001644, XOT00001645, XOT00001646, XOT00001647, XOT00001650, XOT00001655A, XOT00001665, XOT00001671A, XOT00001672A, XOT00001673A, XOT00001676, XOT00001677, XOT00001678, XOT00001680, XOT00001683, XOT00001684, XOT00001685, XOT00001692, XOT00001693, XOT00001694, XOT00001698, XOT00001708, XOT00001710, XOT00001714, XOT00001715, XOT00001722, XOT00001723, XOT00001728A, XOT00001729A, XOT00001736, XOT00001737A, XOT00001750, XOT00001751, XOT00001757, XOT00001758A, XOT00001763, XOT00001764, XOT00001765, XOT00001766, XOT00001767, XOT00001768, XOT00001773, XOT00001841, XOT00001846, XOT00001847, XOT00001848, XOT00001850, XOT00001873A, XOT00001874A, XOT00001883, XOT00001888, XOT00001892, XOT00001893, XOT00001909A, XOT00001925, XOT00001984, XOT00002009, XOT00002020, XOT00002085A, XOT00002136, XOT0012, XOT0015, XOT0016, XOT0017, XOT0018, XOT0019, XOT0020, XOT0021, XOT0022, XOT0023, XOT0032, XOT0033, XOT0034, XOT0035, XOT0036, XOT0038, XOT0042, XOT0044, XOT0052, XOT0053, XOT0054, XOT0055, XOT0056, XOT1005, XOT1006, XOT1007, XOT1008, XOT1009, XOT1012, XOT1027, XOT1032, XOT1034, XOT1040, XOT1049, XOT1051A, XOT1052A, XOT1056, XOT1057, XOT1063, XOT1064, XOT1067, XOT1068, XOT1069, XOT1076, XOT1079, XOT1085, XOT1086, XOT1087, XOT1105, XOT1108, XOT1111, XOT1112, XOT1120, XOT1121, XOT1122, XOT1123, XOT1140, XOT1147, XOT1153, XOT1154, XOT1155, XOT1156, XOT1157, XOT1162, XOT1163, XOT1164, XOT1165, XOT1166, XOT1167, XOT1177, XOT1178, XOT1184, XOT1188, XOT1189, XOT1190, XOT1191, XOT1193, XOT1194, XOT1195, XOT1196, XOT1197, XOT1210, XOT1211, XOT1212, XOT1224, XOT1225, XOT1226, XOT1231, XOT1241, XOT1244, XOT1245, XOT1246, XOT1251, XOT1255, XOT1256, XOT1257, XOT1271, XOT1275, XOT1276, XOT1277, XOT1278, XOT1279, XOT1280, XOT1290, XOT1291, XOT1292, XOT1296, XOT1304, XOT1305, XOT1318, XOT1321, XOT1324, XOT1333, XOT1334, XOT1335, XOT1336, XOT1340, XOT1347, XOT1352, XOT1353, XOT1356, XOT1357, XOT1369, XOT1372, XOT1373, XOT1374, XOT1375, XOT1376, XOT1384, XOT1395, XOT1397, XOT1399, XOT1401, XOT1402, XOT1403, XOT1406, XOT1415, XOT1416, XOT1417, XOT1418, XOT1419, XOT1423, XOT1426, XOT1427, XOT1436, XOT1437, XOT1438, XOT1439, XOT1440, XOT1441, XOT1442, XOT1447, XOT1449, XOT1450, XOT1451, XOT1456, XOT1457, XOT1458, XOT1459, XOT1463, XOT1469, XOT1475A, XOT1477, XOT1483, XOT1485A, XOT1489, XOT1490, XOT1491, XOT1492, XOT1499, XOT1501, XOT1502, XOT1503, XOT1505, XOT1506, XOT1507, XOT1508, XOT1514A, XOT1520A, XOT1527, XOT1531, XOT1532, XOT1543, XOT1544, XOT1547, XOT1556, XOT1557, XOT1560A, XOT1565.
Intended Use of this Product: Indicated for use in the removal of emboli and thrombi during embolectomy and/or thrombectomy. It can also be used for catheter placement over a guidewire, vessel occlusion, fluid infusion and/or aspiration.		
Pruitt® Occlusion Catheter	2103-36, 2103-46, 2103-56, E2103-36, E2103-46, E2103-56	POC1830, POC1840, POC1843, POC1846, POC1847, POC1850, POC1851, POC1856, POC1857, POC1858, POC1860, POC1864, POC1869, POC1873, POC1874, POC1879, POC1890, POC1891, POC1893, POC1897, POC1899, POC1904, POC1912, POC1925, POC1926, POC1927, POC1928, POC1931, POC1932, POC1933, POC1934, POC1935, POC1937, POC1938, POC1939, POC1940, POC1941, POC1942, POC1945, POC1958, POC1962, POC1968, POC1969, POC1971, POC1972, POC1973, POC1974, POC1976, POC1977, POC1979, POC1980, POC1983, POC1997, POC2002, QPO1001, QPO1003, QPO1017, QPO1020, QPO1022, QPO1023, QPO1039, QPO1040, QPO1044, QPO1046, QPO1049, QPO1053, QPO1055, QPO1057, QPO1060, QPO1064, QPO1065, QPO1066, QPO1069, QPO1076, QPO1077, QPO1078, QPO1090, QPO1099, QPO1100, QPO1112, QPO1113, QPO1121, QPO1123, QPO1127, QPO1128, QPO1132, QPO1133, QPO1138, QPO1139, QPO1140, QPO1144, QPO1151, QPO1152, QPO1156, QPO1164, QPO1167, XPO00001036, XPO00001038, XPO00001045, XPO00001051, XPO00001053, XPO00001054, XPO00001061, XPO00001066, XPO00001070, XPO00001077, XPO00001081, XPO00001095, XPO00001098, XPO1001, XPO1002, XPO1004, XPO1011, XPO1013, XPO1018.




Account number: 6213735

CAPA 2025-007

Product Line	Catalog Numbers	Lot Numbers
Intended Use of this Product: <i>The occlusion of vessels both arterial and venous for the control of bleeding.</i>		
Pruitt® Irrigation Occlusion Catheter	2102-09, E2102-09	PIO1454, PIO1455, PIO1457, PIO1458, PIO1463, PIO1466, QPI1001, QPI1006, QPI1012, QPI1013, QPI1014, QPI1015, QPI1016, QPI1018, QPI1025, QPI1028, QPI1029, QPI1030, QPI1031, QPI1032, QPI1037, QPI1038, QPI1040, QPI1042, QPI1043, QPI1045, QPI1049, QPI1050, QPI1053, QPI1054, QPI1055, QPI1060, QPI1062, QPI1064, XPI00001121, XPI00001123, XPI00001129, XPI00001130, XPI00001131, XPI00001141, XPI00001142, XPI00001145, XPI00001161, XPI00001163, XPI1001, XPI1069, XPI1070, XPI1071, XPI1074, XPI1075, XPI1077, XPI1078, XPI1083, XPI1087, XPI1090, XPI1091, XPI1094, XPI1095, XPI1100, XPI1101, XPI1102, XPI1103, XPI1106, XPI1107, XPI1108, XPI1109.
Intended Use of this Product: <i>To temporarily occlude vessels for the control of bleeding. To access the vessel lumen distal to the point of occlusion.</i>		

(3) Reason for the Voluntary Recall

During internal product testing, we observed some packages to have incomplete seals (sterile barrier) as depicted in the picture below. Seal integrity testing has found a compromised package on approximately 1% of devices, which could allow microorganisms to enter and contaminate the device, leading to potential infection transmitted by the compromised device.

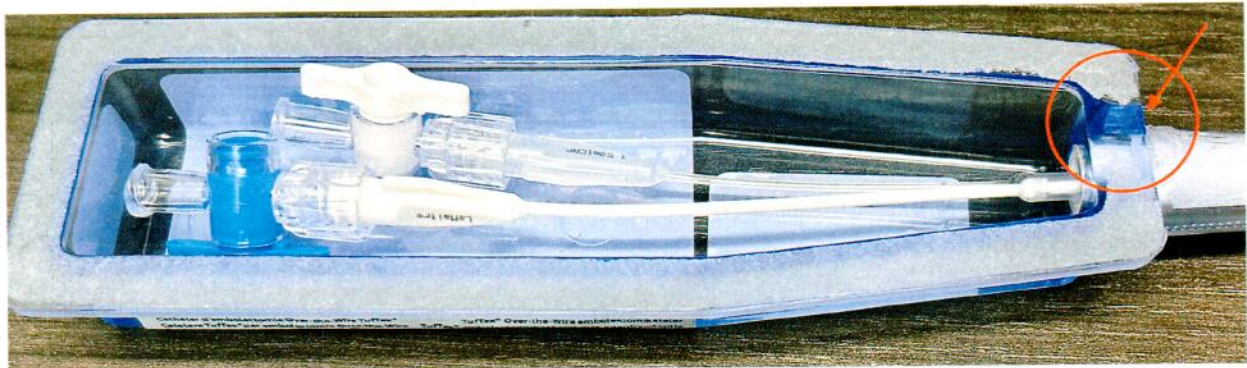

(4) Risk to Health:

If the sterile barrier of the packaged device is broken, there is a risk of contamination of the device and subsequent risk of infection if used. If identified prior to use, there is a risk of delayed treatment due to locating another usable unit. So far, there have been no complaints about this issue and no reports of infection.

(5) Actions to be taken by the Customer/User:

- Do not use the devices. Check your inventory against the list of lots in this letter. Immediately quarantine any recalled devices.
- Complete the form at the end of this letter. **Please note that you must return the form even if you have no devices in inventory.**
- Scan the reply form and send it to recalls@lemaitre.com. (If you have checked your inventory and have none of the recalled devices, you may simply email us with that information to recalls@lemaitre.com.)
 - Alternatively, use the following link to register and complete the reply form electronically:

<https://recall.trievr.com/app/#/?actn=req&k1=10334&>

	Account number: 6213735 CAPA 2025-007
---	--

- Upon receipt of your completed form, a LeMaitre Vascular representative will contact you to provide return instructions.
- When a recalled device has been returned to LeMaitre, a customer service representative will generate a replacement order for the received device(s).
- If you have transferred devices to another facility, please forward a copy of this letter to them.

(6) Product and Distribution Information:

Refer to the tables mentioned earlier in the letter.

(7) Type of Action by the Company:

LeMaitre Vascular has opened a corrective action to address the root causes of the issue. Some actions will be short term, and others will require extensive validation and will therefore be longer term measures.

(8) OTHER INFORMATION:

Authorized by:

Jeffrey Oddy
Sr. Engineer, Post Market Surveillance
Recall Coordinator
Joddy@lemaitre.com

If you have any questions, contact **Jeffrey Oddy** Monday through Friday, 9:00 AM to 5:00 PM, ET.

Adverse reactions or quality problems experienced with the use of this product may be reported to the FDA's MedWatch Adverse Event Reporting program either online, by regular mail or by fax.



Account number: 6213735

CAPA 2025-007

24 Apr 2025 TufTex OTW / POC / PIOC Recall

***This form must be returned even if you have zero devices in inventory.
Email completed form to recalls@lemaitre.com.***

Account #*	Customer Name*	Address
6213735	MOUNT SINAI SOUTH NASSAU	2445 OCEANSIDE ROAD OCEANSIDE, NY 11572

**If you are not the customer listed here, please list your facility information.*

Contact Name	Contact Email	Contact Phone
Signature and Date:		

Contact Information for Ongoing and Future Communications. Input your hospital's risk management contact information below (e.g., riskmanagement@xyzhospital.org). Please do not include an email address for an individual.

Contact Name	Contact Email	Contact Phone

I have read and understand the recall instructions provided in this letter. Yes ☐ No ☐

Any adverse events associated with recalled product(s)? Yes ☐ No ☐

If yes, please explain:

Do you have any recalled devices at your facility? ☐ Yes ☐ No

If Yes, please complete the table below.

- If you have checked your inventory and have no recalled devices, you may simply email recalls@lemaitre.com to indicate that "I have checked our inventory at 6213735, MOUNT SINAI SOUTH NASSAU and we have none of the recalled devices."

NOTE: Distributors must complete the entire form.

- Alternatively, use the following link to register and complete the reply form electronically:
<https://recall.trievr.com/app/#/?actn=reg&k1=10334&>

- **If you have transferred affected devices to another facility, please send them a copy of this recall letter. If possible: list the facility information, including contact information, including contact information.**

REF (Catalog) #	LOT #	QUANTITY ON HAND

	Account number: 6213735
	CAPA 2025-007

REF (Catalog) #	LOT #	QUANTITY ON HAND

ADDRESS TO WHICH REPLACEMENT DEVICES SHOULD BE SENT :

FOR DISTRIBUTORS ONLY:

- ☐ I have checked my stock and have quarantined inventory consisting of ____ units.
- ☐ I identified and notified all of my customers that are affected by this recall.
- ☐ If the product was distributed outside the US, I have notified that country's medical device regulatory agency about this recall.
- ☐ I did not notify the regulatory agency. The rationale is listed below.

Rationale:

Name / Title:	
Telephone:	
Email Address:	

Please scan the completed form and email it to recalls@lemaitre.com.

Thank you for your cooperation.

This section is for LeMaitre use only:

RMA #		REPLACEMENT ORDER #	
--------------	--	----------------------------	--