

A303414A

Correct



Incorrect



A942216

Correct



Incorrect



Clinical Impact:

The current issue may result from a user thinking they are using a silicone catheter, when in fact the tray contains a natural rubber latex catheter.

-or-

The current issue may result from a user thinking they are using a tray containing a natural rubber latex catheter with an infection control coating, when in fact the tray contains a silicone hydrogel coated catheter with antimicrobial coating.

-or-

The current issue may result from a user thinking they are using a tray containing a natural rubber latex catheter with an infection control coating, when in fact the tray contains a silicone hydrogel coated catheter without antimicrobial coating.


In the event a mislabeled tray is used, risk to patients may stem from varying degrees of harm from minor injury to injury that requires medical intervention to treat, such as infection and/or local and/or systemic allergic reaction (anaphylaxis).

Distribution Timeframe and Label Photo:

The affected product was distributed from: 28Dec2024 to 14Mar2025

Sample Label Photos:


ProductLabel

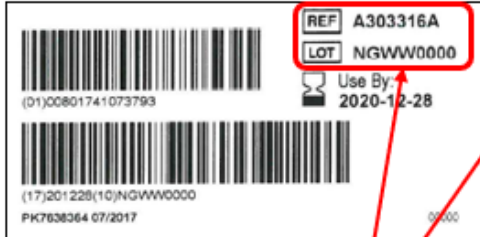

Case/Carton Label



Location of REF# and LOT#

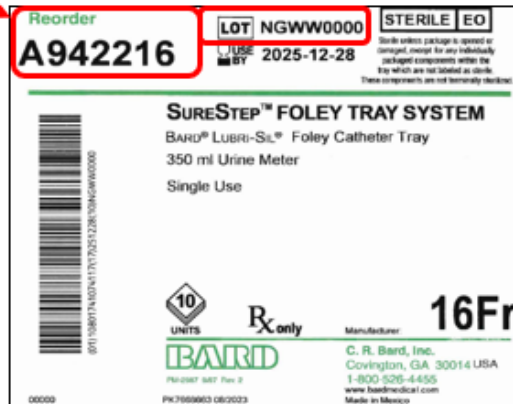

Location of REF# and LOT#



Please Take the Following Actions:

1. **Immediately Discontinue Distribution** of all affected product codes, immediately complete the destruction and disposal of remaining inventory in your possession in accordance with your local facility's process.
2. Identify all customers within your distribution network that purchased any affected product as defined in this notification. Provide a copy of the attached customer letter to all customers to advise them of this field action notification on BD's behalf. Attachment (2) (Customer Letter)
3. Share this notice with any users of the product within your distribution network to ensure they are also aware of this Medical Device Recall (Correction).
4. Complete the attached Distributor Response Form and return to the BD contact noted on the form confirming acknowledgement of the recall notification, **whether or not you have any of the impacted product**, so that BD may acknowledge your receipt of this notification.
5. Indicate on the response form the quantity of product identified at your facility and confirm that this inventory was destroyed. Attachment (1) (Distributor Response Form)
6. Please contact your BD representative if you require assistance with credits to be issued for discard-ed products or for replacement product options. You may also contact the number below for further assistance.
7. Report any adverse health consequences experienced with the use of this product to BD. Events may also be reported to the FDA's MedWatch Adverse Event Reporting program via:

Web: MedWatch website at www.fda.gov/medwatch

Phone: 1-800-FDA-1088 (1-800-332-1088)

Mail: MedWatch, HF-2, FDA, 5600 Fisher's Lane, Rockville, MD 20852-9787

8. Report any complaints experienced with the use of this product to BD via the North American Regional Complaint Center:
 Phone: 1-844-8BD-LIFE (1-844-823-5433) say "Product Complaints" when prompted
 Mon–Fri 8:00am and 5:00pm CT
 Email: productcomplaints@bd.com

Affected Product:

Product Name	Catalog No.s	Lot Numbers	UDI-DI	Expiration Date DD-MMM-YYYY	Product Package Size
SureStep™ Foley Tray System Bardex® I.C. Complete Care® Infection Control Foley Catheter Tray	A303316A	NGJX6399	Eaches (01)00801741073793 Case (01)10801741073790	31May2027	EA Case
SureStep™ Foley Tray System Lubri-Sil® I.C. Complete Care® Infection Control Foley Catheter Tray	A303414A	NGJY3685	Eaches (01)00801741073816 Case (01)10801741073813	31May2027	EA Case

SureStep™ Foley Tray System BARD® Lubri-Sil® Foley Catheter Tray	A942216	NGJX4320	Eaches (01)00801741074110 Case (01)10801741074117	31May2027	EA Case
--	---------	----------	--	-----------	------------

Actions Taken by BD:

BD is investigating the cause of the issue and will implement corrective action based on the investigation results to prevent re-occurrence.

Contact Information: If you require further assistance please contact:

BD Contact	Contact Information	Areas of Support
BD Medical Information Services	Available 24/7 Select #2 when prompted 1-800-555-7422 or Email: medical.information@bd.com	Product Use Inquiries
BD Technical Support	1-844-8BD-LIFE (1-844-823-5433) Say "Product Complaint" when prompted M-F 8am - 5pm CT or Email: productcomplaints@bd.com	Product Complaints, Technical Support
Recall Customer Response Forms	Email: BDRC38@bd.com Fax No.: 312-949-0574	Mail for receipt of customer response forms

BD is committed to advancing the world of health. Our primary objectives are patient and user safety and providing you with quality products. We apologize for any inconvenience this issue may have caused and thank you in advance for helping us to resolve this matter as quickly and effectively as possible.

Sincerely,



Gerald Denny
Vice President, Medical Affairs
BD UCC



Elizabeth Gaipa
Vice President, Quality Management
BD UCC

Enclosures:

Attachment 1 – Distributor Response Form
Attachment 2 – Customer Letter

**DISTRIBUTOR RESPONSE FORM
(UCC-25-5282, UCC-25-5319)**

SureStep™ Foley Tray Systems: Bardex® I.C. Complete Care® Infection Control Foley Catheter Tray; Lubri-Sil® I.C. Complete Care® Infection Control Foley Catheter Tray; BARD® Lubri-Sil® Foley Catheter Tray

Please assist BD by acknowledging this field action using one of the following methods:

Email: BDRC38@bd.com

Fax No.: 312-949-0574

Facility: _____
Please use full, current facility name. Do not use initials.

Street Address: _____

City: _____ State: _____ Zip: _____

Field action Response Form Completed By:	
Name:	
Title:	
Telephone No.	
Fax No.	
Email Address	

Please check all that apply:

- ☐ I have read and understood the attached notice and taken appropriate actions.
- ☐ We do not have any of the affected product(s) on hand.
- ☐ We have destroyed all affected product indicated below as available inventory at the time of receipt of this notification and request credit for the following units.

Product Name	Catalog No.	Lot No.	No. of Units (Eaches/Cases)
SureStep™ Foley Tray System Bardex® I.C. Complete Care® Infection Control Foley Catheter Tray	A303316A	NGJX6399	
SureStep™ Foley Tray System Lubri-Sil® I.C. Complete Care® Infection Control Foley Catheter Tray	A303414A	NGJY3685	
SureStep™ Foley Tray System BARD® Lubri-Sil® Foley Catheter Tray	A942216	NGJX4320	

Please assist BD with assuring these communications are delivered to the appropriate person/function within your facility if that is not you.

Person/function responsible for the receipt and management of all recalls/corrections at your facility:

Name: _____

Phone: _____

Email: _____

Fax: _____

Attachment 1

URGENT Medical Device Recall

SureStep™ Foley Tray Systems: Bardex® I.C. Complete Care® Infection Control Foley Catheter Tray; Lubri-Sil® I.C. Complete Care® Infection Control Foley Catheter Tray; BARD® Lubri-Sil® Foley Catheter Tray

July 17, 2025

For the Attention of: Recall Coordinator, Director of Nursing, Director of Purchasing, Director of Risk Management

The purpose of this letter is to inform you that C.R. Bard® Urology and Critical Care, a wholly owned subsidiary of Becton, Dickinson and Company (BD) is voluntarily initiating a correction for the devices listed below. All products listed in the Affected Product section are impacted by this recall and includes the following products:

- **SureStep™ Foley Tray System Bardex® I.C. Complete Care® Infection Control Foley Catheter Tray**
- **SureStep™ Foley Tray System Lubri-Sil® I.C. Complete Care® Infection Control Foley Catheter Tray**
- **SureStep™ Foley Tray System BARD® Lubri-Sil® Foley Catheter Tray**

Description of the problem:

The incorrect product information insert sheet was included with the product which does not indicate the correct information including the material used to produce it, Silicone or Natural Rubber Latex, and presence of related antimicrobial coated components. See photos below.

A303316A

Correct



Incorrect



A303414A
Correct



Incorrect



A942216
Correct



Incorrect



Clinical Impact:

The current issue may result from a user thinking they are using a silicone catheter, when in fact the tray contains a natural rubber latex catheter.

-or-

The current issue may result from a user thinking they are using a tray containing a natural rubber latex catheter with an infection control coating, when in fact the tray contains a silicone hydrogel coated catheter with antimicrobial coating.

-or-

The current issue may result from a user thinking they are using a tray containing a natural rubber latex catheter with an infection control coating, when in fact the tray contains a silicone hydrogel coated catheter without antimicrobial coating.

In the event a mislabeled tray is used, risk to patients may stem from varying degrees of harm from minor injury to injury that requires medical intervention to treat, such as infection and/or local and/or systemic allergic reaction (anaphylaxis).

Advice for Clinical Users:


It is recommended that clinicians discontinue the use of impacted devices and find alternative catheters with the appropriate characteristics required for their patients to continue therapy and avoid potential health risks related to this issue.

Distribution Timeframe and Sample Label Photo:

The affected product was distributed from: 17Jan2025 to 12Mar2025

Sample Label Photos:

ProductLabel

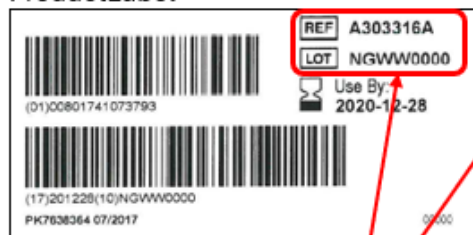

Case/Carton Label



Location of REF# and LOT#



Location of REF# and LOT#



Please Take the Following Actions:

1. **Immediately Discontinue Use.**
2. Please check all inventory locations within your institution for affected BD (C.R. Bard, Inc.) product listed in the Affected Product section below.
3. Immediately quarantine and discard all devices within your facility's control per your facility's procedures.
4. Share this notice with any users of the product within your facilities to ensure they are also aware of this Medical Device Recall (Correction).
5. If you have purchased this product through a distributor, please contact this distributor to request appropriate credit or order alternate replacement product.
6. Complete and return the attached Customer Response Form to the BD contact noted on the form confirming acknowledgement of the correction notification, **whether or not you have any affected product**, so that BD may confirm your receipt of this notification and provide credit or alternate product. Attachment (1) (Customer Response Form)
7. Report any adverse health consequences experienced with the use of this product to BD. Events may also be reported to the FDA's MedWatch Adverse Event Reporting program via:

Web: MedWatch website at www.fda.gov/medwatch

Phone: 1-800-FDA-1088 (1-800-332-1088)

Mail: MedWatch, HF-2, FDA, 5600 Fisher's Lane, Rockville, MD 20852-9787

8. Report any complaints experienced with the use of this product to BD via the North American Regional Complaint Center:
Phone: 1-844-8BD-LIFE (1-844-823-5433) say "Product Complaints" when prompted Mon–Fri 8:00am and 5:00pm CT
Email: productcomplaints@bd.com

Affected Product:

Product Name	Catalog No.s	Lot Numbers	UDI-DI	Expiration Date DD-MMM-YYYY	Product Package Size
SureStep™ Foley Tray System Bardex® I.C. Complete Care® Infection Control Foley Catheter Tray	A303316A	NGJX6399	Eaches (01)00801741073793 Case (01)10801741073790	31May2027	EA Case
SureStep™ Foley Tray System Lubri-Sil® I.C. Complete Care® Infection Control Foley Catheter Tray	A303414A	NGJY3685	Eaches (01)00801741073816 Case (01)10801741073813	31May2027	EA Case

SureStep™ Foley Tray System BARD® Lubri-Sil® Foley Catheter Tray	A942216	NGJX4320	Eaches (01)00801741074110 Case (01)10801741074117	31May2027	EA Case
--	---------	----------	--	-----------	------------

Actions Taken by BD:

BD is investigating the cause of the issue and will implement corrective action based on the investigation results to prevent re-occurrence.

Contact Information: If you require further assistance please contact:

BD Contact	Contact Information	Areas of Support
BD Medical Information Services	Available 24/7 Select #2 when prompted 1-800-555-7422 or Email: medical.information@bd.com	Product Use Inquiries
BD Technical Support	1-844-8BD-LIFE (1-844-823-5433) Say "Product Complaint" when prompted M-F 8am - 5pm CT or Email: productcomplaints@bd.com	Product Complaints, Technical Support
Recall Customer Response Forms	Email: BDR38@bd.com Fax No.: 312-949-0574	Mail for receipt of customer response forms

BD is committed to advancing the world of health. Our primary objectives are patient and user safety and providing you with quality products. We apologize for any inconvenience this issue may have caused and thank you in advance for helping us to resolve this matter as quickly and effectively as possible.

Sincerely,



Gerald Denny
Vice President, Medical Affairs
BD UCC



Elizabeth Gaipa
Vice President, Quality Management
BD UCC

Enclosures:

Attachment 1 – Customer Response Form

**CUSTOMER RESPONSE FORM
(UCC-25-5282, UCC-25-5319)**

**SureStep™ Foley Tray Systems: Bardex® I.C. Complete Care® Infection Control
Foley Catheter Tray; Lubri-Sil® I.C. Complete Care® Infection Control Foley
Catheter Tray; BARD® Lubri-Sil® Foley Catheter Tray**

Please assist BD by acknowledging this field action using one of the following methods:

Email: BDRC38@bd.com

Fax No.: 312-949-0574

Facility: _____
Please use full, current facility name. Do not use initials.

Street Address: _____

City: _____ State: _____ Zip: _____

Field action Response Form Completed By:	
Name:	
Title:	
Telephone No.	
Fax No.	
Email Address	

Please check all that apply:

- ☐ I have read and understood the attached notice and taken appropriate actions.
- ☐ We do not have any of the affected product(s) on hand.
- ☐ I have destroyed all affected product indicated below as available inventory at the time of receipt of this notification and request credit for the following units.
- ☐ I have destroyed all affected product indicated below as available inventory at the time of receipt of this notification and request replacement for the following units.

Product Name	Catalog No.	Lot No.	No. of Units (Eaches/Cases)
SureStep™ Foley Tray System Bardex® I.C. Complete Care® Infection Control Foley Catheter Tray	A303316A	NGJX6399	
SureStep™ Foley Tray System Lubri-Sil® I.C. Complete Care® Infection Control Foley Catheter Tray	A303414A	NGJY3685	
SureStep™ Foley Tray System BARD® Lubri-Sil® Foley Catheter Tray	A942216	NGJX4320	

Please assist BD with assuring these communications are delivered to the appropriate person/function within your facility if that is not you.

Person/function responsible for the receipt and management of all recalls/corrections at your facility:

Name: _____

Phone: _____

Email: _____

Fax: _____

Attachment 1